

Multiple High-Resolution Serum Proteomic Features for Ovarian Cancer Detection

Background

[1001] Serum proteomic pattern analysis by mass spectrometry (MS) is an emerging technology that is being used to identify biomarker disease profiles. Using this MS-based approach, the mass spectra generated from a training set of serum samples is analyzed by a bioinformatic algorithm to identify diagnostic signature patterns comprised of a subset of key mass-to-charge (m/z) species and their relative intensities. Mass spectra from unknown samples are subsequently classified by likeness to the pattern found in mass spectra used in the training set. The number of key m/z species whose combined relative intensities define the pattern represent a very small subset of the entire number of species present in any given serum mass spectrum.

[1002] The feasibility of using MS proteomic pattern analysis for the diagnosis of ovarian, breast, and prostate cancer has been demonstrated. While investigators have used a variety of different bioinformatic algorithms for pattern discovery, the most common analytical platform is comprised of a low-resolution time-of-flight (TOF) mass spectrometer where samples are ionized by surface enhanced laser desorption/ionization (SELDI), a ProteinChip array-based chromatographic retention technology that allows for direct mass spectrometric analysis of analytes retained on the array.

[1003] Ovarian cancer is the leading cause of gynecological malignancy and is the fifth most common cause of cancer-related death in women. The American Cancer Society estimates that there will be 23,300 new cases of ovarian cancer and 13,900 deaths in 2002. Unfortunately, almost 80% of women with common epithelial ovarian cancer are not diagnosed until the disease is advanced in stage, i.e., has spread to the upper abdomen (stage III) or beyond (stage IV). The 5-year survival rate for these women is only 15 to 20%, whereas the 5-year survival rate for ovarian cancer at stage I approaches 95% with surgical intervention. The early diagnosis of ovarian cancer, therefore, could dramatically decrease the number of deaths from this cancer.

[1004] The most widely used diagnostic biomarker for ovarian cancer is Cancer Antigen 125 (CA 125) as detected by the monoclonal antibody OC 125. Though 80% of patients with ovarian cancer possess elevated levels of CA 125, it is elevated in only 50-60% of patients at stage I, lending it a positive-predictive value of 10%. Moreover, CA 125 can be elevated in other non-gynecologic and benign conditions. A combined strategy of CA 125 determination with ultrasonography increases the positive-predictive value to approximately 20%.

[1005] Low molecular weight serum proteomic patterns from low-resolution SELDI-TOF MS data can distinguish neoplastic from non-neoplastic disease within the ovary. See Petticon, R. P. III, *et al.* Use of proteomic patterns in serum to identify ovarian cancer. *The Lancet* 359, 572-577 (2002). The proteomic patterns can be identified by application of an artificial intelligence bioinformatics tool that employs an unsupervised system (self-organizing cluster mapping) as a fitness test for a supervised system (genetic algorithm). A training set comprised of SELDI-TOF mass spectra from serum derived from either unaffected women or women with ovarian cancer is employed so that the most fit combination of m/z features (along with their relative intensities) plotted in n -space can reliably distinguish the cohorts used in training. The "trained" algorithm is applied to a marked set of samples that resulted in a sensitivity of 100% and a specificity of 93%. This technique is described in more detail in WO 02/06839A2 "A Process for Artificially Directed Molecular Sorting Based on Hidden Patterns from Biological Data" ("Hidden Patterns") the disclosure of which is hereby expressly incorporated herein by reference.

[1006] Although this technique works well, the low-resolution mass spectrometric instrumentation and thus the data that comes from the instrument may limit the attainable reproducibility, sensitivity, and specificity for proteomic pattern analyses for routine clinical use.

Summary

[1007] The protein pattern analysis concept of Hidden Patterns is extended to a high-resolution MS platform to generate diagnostic models possessing higher sensitivities and

specificities on a format that generates more stable spectra, has a true time-of-flight mass accuracy, and is inherently more reproducible, machine-to-machine, and day-to-day because of the increase in mass accuracy. Sera from a large, well-controlled ovarian cancer screening trial were used and proteomic pattern analysis was conducted on the same samples on two mass spectral platforms differing in their effective resolution and mass accuracy. The data was analyzed to also rank the sensitivity and specificity of the series of diagnostic models that emerged.

[1008] The spectra from a high-resolution and a low-resolution mass spectrometer with the same patient sera samples applied and analyzed on the same SELDI ProteinChip arrays were compared. Although the higher resolution mass spectra may generate more distinguishable sets of diagnostic features, the increased complexity and dimensionality of data may reduce the likelihood of fruitful pattern discovery. Diagnostic proteomic feature sets can be discerned within the high-resolution spectra from the clinically relevant spectral study set, and the modeling outcomes between the two instrument platforms can be compared. The number and character of the diagnostic models emerging from data mining operations can be tracked. Serum proteomic pattern analysis can be used for the generation of multiple, highly accurate models using a hybrid quadrupole, time-of-flight (Qq-TOF) MS for an improved early diagnosis of ovarian cancer.

Brief Description of the Figures

[1009] FIGS. 1A and 1B compare the mass spectra from control serum prepared on a WCCX2 ProteinChip array and analyzed with a PBS-II-TOF (panel A) or a Qq-TOF (panel B) mass spectrometer.

[1010] FIGS. 2A and 2B show histograms representing the testing results of sensitivity (2A) and specificity (2B) of 108 models for MS data acquired on either a Qq-TOF or a PBS-II-TOF mass spectrometer.

[1011] FIGS. 3A and 3B show histograms representing the testing and blinded validation results of sensitivity (3A) and specificity (3B) of 108 models for MS data acquired on either a Qq-TOF or a PBS-II-TOF mass spectrometer.

[1012] FIGS. 4A and 4B compare SELDI Qq-TOF mass spectra of serum from an unaffected individual (4A) and an ovarian cancer patient (4B).

Detailed Description

Analysis of Serum Samples

[1013] A total of 248 serum samples were provided from the National Ovarian Cancer Early Detection Program (NOEDP) clinic at Northwestern University Hospital (Chicago, Illinois). The samples were processed and their proteomic patterns acquired by MS as described below in the description of the methods used. The serum samples in the present study were analyzed on the same protein chip arrays by both a PBS-II and a Qq-TOF MS fitted with a SELDI ProteinChip array interface. While the spectra acquired from both instruments are qualitatively similar, the higher resolution afforded by the Qq-TOF MS is apparent from FIG. 1. This increased resolution allows species close in m/z unresolved by the PBS-II-TOF MS to be distinctly observed in the Qq-TOF mass spectrum. Indeed, simulations demonstrate the ability of the Qq-TOF MS (routine resolution ~ 4000) to completely resolve species differing in m/z of only 0.273 (e.g., at m/z 3000) whereas complete resolution of species with the PBS-II-TOF MS (routinely resolution ~ 150) is only possible for species that differ by m/z of 20 (simulation not shown).

[1014] The mass spectra were analyzed using the ProteomeQuestTM bioinformatics tool employing ASCII files consisting of m/z and intensity values of either the PBS-II-TOF or the Qq-TOF mass spectra as the input. The mass spectral data acquired using the Qq-TOF MS were binned to precisely define the number of features in each spectrum to 7084 with each feature being comprised of a binned m/z and amplitude value. The algorithm examines the data to find a set of features of precise binned m/z values whose combined, normalized relative intensity values in a space best segregate the data derived

from the training set. Mass spectra acquired on the Q4-TOF and the YBS-II TOF instruments from the same sample sets were restricted to the m/z range from 700 to 11,893 for direct comparison between the two platforms. The entire set of spectra acquired from the serum samples was divided into three sets: a) a training set that is used to discover the hidden diagnostic patterns; b) a testing set; and c) a validation set. With this approach, only the normalized intensities of the key m/z values were identified using the training set were used to classify the testing and validation sets, and the algorithm had not previously "seen" the spectra in the testing and validation sets.

[1015] The training set was comprised of serum from 28 unaffected women and 56 women with ovarian cancer. The training and testing set mass spectrometry analyzed by the blindfolded algorithm: to generate a series of models under the following set modeling parameters: a) a randomly space of 85%, 90%, or 95% likelihood for disease classification; b) a feature set size of 5, 10, or 15 random m/z values whose combined intensities comprise each pattern; and c) a learning rate of 0.1%, 0.2%, or 0.3% for pattern generation by the genetic algorithm. Four sets of randomly generated models for each of the 27 permutations were derived and queried with the same test set. Sensitivity and specificity testing results for each of the 108 models (four rounds of training for each of the 27 permutations) were generated, as shown in FIGS. 2A and 2B. These results demonstrate that the Q4-TOF MS did produce better results than the lower resolution version ($P < 0.00001$, using the exact Cochran-Armitage test (GSA) great-Armitage test, Data Analysis New York: John Wiley and Sons (1999)) for serum throughout a range of modeling conditions.

[1016] The ability to generate the best performing models for testing and validation was statistically evaluated as multiple models were generated and ranked using the entire range of the modeling parameters above. Models from the training set were validated using a testing set consisting of 31 unaffected and 63 ovarian cancer serum samples. To further validate the ability to diagnose ovarian cancer, a set of blinded sample mass spectra consisting of an additional 37 normal and 40 ovarian cancer serum mass spectra were tested against the model found in training previously discussed. As shown in FIGS. 3A and 3B, the results show the ability of the mass spectra from the higher

resolution Q4-TOF MS to generate statistically significant ($P < 0.00001$) superior models over the lower resolution YBS-II mass spectra.

[1017] Fifteen models were found that were 100% sensitive in their ability to correctly discriminate unaffected women from those suffering from ovarian cancer, that were 100% specific to discriminating women in the test set, and at least 97% specific in the validation set. These models are shown in Appendix A, and identified as Model 1 through Model 15. Of these models, four were found that were both 100% sensitive and specific for both sets (Models 4, 9, 10, and 15).

[1018] Appendix A identifies for each model the following information: First, the specificity and sensitivity for each model is shown for the Test set and for the Validity set. The number of samples for which the model correctly grouped women with a "Normal State" (i.e. not having ovarian cancer) and with an "Ovarian Cancer State" is then shown for each of the test and validity sets, compared to the total number of samples in the corresponding sets. For example, in Model 1, the model correctly identified 36 of the 37 women as having a normal state in the Validity set.

[1019] Finally, for each model a table is set forth showing the constituent "patterns" comprising the model. Each pattern corresponds to a point, or node, in the N-dimensional space defined by the N m/z values (or "features") included in the model. Thus, each pattern is a set of features, each feature having an amplitude. Appendix A therefore shows for each model a table containing the constituent patterns, each pattern being in a row identified by a "Node" number. The table also includes columns for the constituent features of the patterns, with the m/z value for each pattern identified at the top of the column. The amplitudes are shown for each feature, for each pattern, and are normalized to 1.0. The remaining four columns in each table are labeled "Count," "State," "Status," and "Error." "Count" is the number of samples in the Training set that correspond to the identified node. "State" indicates the state of the node, where 1 indicates diseased (in this case, having ovarian cancer) and 0 indicates normal (not having the disease). "Status" is the sum of the state values for all of the correctly classified members of the indicated node, while "Error" is the number of incorrectly

classified members of the indicated node. Thus, for node 5 in Model 1, 13 samples were assigned to the node, whereas 11 samples were actually diseased. Statistics in lines 11 (number 13) and Error is 2.

[1020] Examination of the key *m/z* features that comprise the four best performing models (Models 4, 9, 10, and 15) reveals certain features (i.e., contained within *m/z* bins 7060-121, 8603-678 and 8706-065) that are consistently present as classifiers in those models.

[1021] Although the proteomic patterns generated from both healthy and cancer patients using the Q-TOF MS are quite similar (as seen by comparing FIGS. 4A to 4B), careful inspection of the raw mass spectra reveals that peaks within the binned *m/z* values 7060-121 and 8603-678 are differentially abundant. In a selection of the serum samples obtained from ovarian cancer patients as compared to unaffected individuals and that the features that the ProteomicQuest™ software selected are "real" features and not noise. The insets in FIGS. 4A and 4B show expanded *m/z* regions highlighting significant intensity differences of the peaks in the *m/z* bins 7060-121 and 8603-678 (indicated by brackets) identified by the algorithm as belonging to the optimum discriminatory pattern. These results indicate these MS peaks originate from species that may be consistent indicators of the presence of ovarian cancer. The ability to distinguish sera from an unaffected individual or an individual with ovarian cancer based on a single serum proteomic *m/z* feature alone, however, is not possible because the entire serum proteome, while a single key *m/z* species is insufficient to globally distinguish all of the unaffected and ovarian cancer patients taken together, the combined peak intensities of key ions does allow the two data sets to be completely distinguished.

[1022] The four best performing models that are 100% sensitive and specific for the blinded testing and validation tests were chosen for further analysis. Table 1 shows bioinformatic classification results of serum samples from masked testing and validation sets by proteomic pattern classification using the best performing models.

	Actual	Predicted (%)
Benign / Unaffected	66	66 (100)
Ovarian Cancer Stage I	22	22 (100)
Ovarian Cancer Stage II, III, IV	51	51 (100)

Each of these models was able to successfully diagnose the presence of ovarian cancer in all of the serum samples from affected women. Further, no false positive or false negative classifications occurred with these best performing models.

Discussion

[1023] A limitation of individual cancer biomarkers is the lack of sensitivity and specificity when applied to large heterogeneous populations. Biomarker pattern analysis seeks to overcome the limitation of individual biomarkers. Serum proteomic pattern analysis can provide new tools for early diagnosis, therapeutic monitoring and outcome analysis. Its usefulness is enhanced by the ability of a selected set of features to transcend the biologic heterogeneity and methodological background "noise." This diagnostic goal is aided by employing a genetic algorithm coupled with a self-organizing cluster analysis to discover diagnostic subsets of *m/z* features and their relative intensities contained within high-resolution Q-TOF mass spectral data.

[1024] It is believed that diagnostic serum proteomic feature sets exist within constellations of small proteins and peptides. A given signature pattern reflects changes in the physiologic or pathologic state of a target tissue. With regard to cancer markers, it is believed that serum diagnostic patterns are a product of the complex, tumor-host immunoreaction. It is thought, likely, that the set of diagnostic features is partially derived from multiple modified host proteins rather than consisting exclusively from the cancer cells. The biomarker profile may be amplified by tumor-host interactions. This amplification includes, for example, the generation of peptide cleavage products by tumor or host processes. There may exist multiple dependent or independent sets of proteins/peptides that reflect the underlying tissue pathology. Hence, the disease related proteomic pattern information content in blood might be richer than previously anticipated. Rather than a single "best" feature set, multiple proteomic feature sets may exist that achieve highly accurate discrimination and hence diagnostic power. This possibility is supported by the data described above.

(1025) The low molecular weight serum proteinome is an unexplored archive, even though this is the mass region where MS is best suited for analysis. It is thought likely that disease-associated species are comprised of low molecular weight peptides/proteins species that vary in mass by as little as a few Daltons. Thus a higher resolution mass spectrometer would be expected to discriminate and discover patterns not resolvable by a lower resolution instrument. The spectra produced by a Q-TOF MS were compared to that of the Capgemon PBS-II TOF MS. The routine resolution obtained is in excess of 8000 ($m/z = 1500$) for the Q-TOF MS and 150 ($m/z = 1500$) for the PBS-II TOF mass spectrometer. A SEED source was used so that both instruments analyzed the same sample on distinct regions of the protein chip array ball surface. While the overall spectral profile is similar, a single peak on the PBS-II TOF MS is resolved into a multitude of peaks on the Q-TOF MS (seen by comparing FIGS 1A and 1B to FIGS. 4A and 4B). Moreover, the inherent increase in mass accuracy by higher resolution instrumentation that has uncoupled the mass analyzer from the source will provide for cleaner spectra as this will suppress confounding, ineluctable ions, generate spectra with lower mass drift over time and instruments at the same time as generating more complex, highly resolved data.

[1026] In the first phase of comparison, problematic problems from mass spectra derived from the same training set and generated on the high and low-resolution mass spectrometers were scrutinized for high overall sensitivity and specificity over a series of modeling constraints in which problems were generated using three different degrees of modeling constraints for the self-organizing clusters to form, three different sets of feature similarity space for the self-organizing clusters to form, three different sets of feature sizes chosen, and three different mutation rates for a total of 27 modeling permutations. Sensitivity and specificity testing results for each of the 108 models (shown in FIGS. 2A and 2B), produced from four rounds of training the set of the 27 permutations, demonstrate that the Qc-TQF MS generated spectra consistently outperformed the lower resolution TOP-MS spectra ($p < 0.00001$) independent of the modeling criteria used.

Since the spectra from the higher resolution platform generate patterns with a higher level of sensitivity and specificity, those spectra could generate more accurate models with a higher degree of sensitivity and specificity - that is, generate the test

diagnostic models. These results were generated using even more stringent criteria, in that an additional nested validation set was employed after testing to determine overall accuracy. The higher resolution spectra consistently produced significantly more accurate models as seen in both the testing and validation studies (as shown in FIGS. 3A and 3B). The models derived from the Qc-TOP MS were consistently more sensitive and specific ($P < 0.00001$) than those from the PBS-II TOP MS. Four models were generated that attained 100% sensitivity and specificity in both testing and validation. The number of key *m/z* values used in classification in the four best diagnostic models ranged from 5 to 9. Three *m/z* bin values were found in two of these four models and two *m/z* bin values were found in three of the four best models. The distinct peaks present in the recurring *m/z* bins 7060, 121, 8503, 6773 and 8706,065 may be good candidates for low molecular weight components in serum that may be key disease progression indicators.

11029] These data support the existence of multiple highly accurate and distinct proteomic feature sets that can accurately distinguish ovarian cancer. To screen for diseases of relatively low prevalence, such as ovarian cancer, a diagnostic test preferably encodes 90% sensitivity and specificity to minimize false positives, while correctly detecting early stage disease when it is present. As discussed above, four models generated using high-resolution Q1-TOF MS data achieved 100% sensitivity and specificity. In blinded testing and validation studies any one of these models were used to correctly classify 22/23 stage-I ovarian cancer, 8/9 stage-II ovarian cancer, stages II, III, and IV and 68/68 benign disease controls.

[1029] Thus, a clinical test could simultaneously employ several combinations of highly accurate diagnostic patterns arising concomitantly from the same data streams, which, taken together, could achieve an even higher degree of accuracy in a screening setting where a diagnostic test will face large population heterogeneity and potential variability in sample quality and handling. Hence, a high-resolution system, such as the Qc-TOR MS employed in this study, is preferred based on the present results,

Methods

[1030] **Serum Samples:** Serum samples were obtained from the National Ovarian Cancer Early Detection Program (NOEDP) clinic at Northwestern University Hospital (Chicago, Illinois). Two hundred and forty eight samples were prepared using a Biomek (Chicago, Illinois). Two hundred and forty eight samples were prepared using a Biomek 2000 robotic liquid handler (Beckman Coulter, Inc., Palo Alto, California). All analyses were performed using ProteinChip weak cation exchange interaction chips (WCX2, Cytogen Biosystems Inc., Fremont, California). A control sample was randomly applied to one spot on each protein array as a quality control for sample preparation and mass spectrometer function. The control sample, SRM 1931A, which is comprised of pooled human sera, was provided by the National Institute of Standards and Technology (NIST).

[1031] **Sample Preparation:** WCX2 ProteinChip arrays were processed in parallel using a Biomek Laboratory workstation (Beckman-Coulter) modified to make use of a ProteinChip array bioprocessor (Cytogen Biosystems Inc.). The bioprocessor holds 12 ProteinChips, each having 8 chromatographic "spots" allowing 96 samples to be processed in parallel. One hundred μ l of 10 mM HCl was applied to the WCX2 protein arrays and allowed to incubate for 5 minutes. The HCl was aspirated, discarded and 100 μ l of distilled, deionized water (ddH₂O) was applied and allowed to incubate for 1 minute. The ddH₂O was aspirated, discarded, and replaced for another minute. One hundred μ l of 10 mM NH₄HCO₃ with 0.1% Triton X-100 was applied to the surface and allowed to incubate for 5 minutes after which the solution was aspirated and discarded. A second application of 100 μ l of 10 mM NH₄HCO₃ with 0.1% Triton X-100 was applied and allowed to incubate for 5 minutes after which the ProteinChip array had residues were aspirated. Five μ l of raw, undiluted serum was applied to each ProteinChip WCX2 that surface and allowed to incubate for 55 minutes. Bead ProteinChip array was washed 3 times with Dulbecco's phosphate buffered saline (PBS) and ddH₂O. For each wash, 150 μ l of either PBS or ddH₂O was sequentially dispensed, mixed by aspirating, and dispensed for a total of 10 times in the bioprocessor after which the solution was aspirated to waste. This wash process was repeated for a total of 6 washes per ProteinChip array half surface. The ProteinChip array half surfaces were vacuum dried to prevent cross contamination when the bioprocessor gases was removed. After removing

the bioprocessor gases, 1.0 μ l of a saturated solution of α -cyano-5-hydroxymethyl acid in 50% (v/v) acetonitrile, 0.5% (v/v) trifluoroacetic acid was applied to each spot on the ProteinChip array twice, allowing the solution to dry between applications.

[1032] **PBS-II Analysis:** ProteinChip arrays were placed in the Protein Biological System II time-of-flight mass spectrometer (PBS-II, Cytogen Biosystems Inc.) and mass spectra were recorded using the following settings: 195 laser shot/spectrum collected in positive mode, laser intensity 220, detector sensitivity 5, detector voltage 1850, and a mass focus of 6000 Da. The PBS-II was externally calibrated using the "All-In-One" peptide mass standard (Cytogen Biosystems, Inc.).

[1033] **Q4-TOF MS Analysis:** ProteinChip arrays were analyzed using a hybrid quadrupole time-of-flight mass spectrometer (QSTAR pulsar I, Applied Biosystems Inc., Framingham, Massachusetts) fitted with a ProteinChip array interface (Cytogen Biosystems Inc., Fremont, California). Samples were ionized with a 337 nm pulsed nitrogen laser (Thornelaser Sciences model VSL-337-ND-S, Walham, Massachusetts) operating at 30 Hz. Approximately 20 mTorr of nitrogen gas was used for collisional ion cooling. Each spectrum represents 100 multi-channel averaged scans (1.667 min acquisition/spectrum). The mass spectrometer was externally calibrated using a mixture of known peptides.

[1034] **Peptide, Protein, and Gene Analysis:** Peptide, protein, and gene analysis was performed by exporting the raw data file generated from the Q4-TOF mass spectrometer into a tab delimited format that generated approximately 350,000 data points per spectrum. The data files were binned using a function of 400 parts per million (ppm) such that all data files possess identical m/z values (e.g., the m/z bin sizes linearly increased from 0.28 at m/z 700 to 4.75 at m/z 12,000). The intensities in each 400 ppm bin were summed. This binning process condenses the number of data points to exactly 7,084 points per sample. The binned spectral data were separated into approximately three equal groups for training, testing and blind validation. The training set consisted of 28 normal and 56 ovarian cancer samples. The models were built on the training set using ProteinQuest™ (Coriolis Systems Inc., Bethesda, Maryland) and validated using the

testing samples, which consisted of 30 normal and 57 ovarian cancer samples. The model was validated using blinded samples, which consisted of 37 normal and 40 ovarian cancer samples. These m/z values that were found to be classifiers were to distinguish between a patient with ovarian cancer from that of an unaffected individual are based on the blinded data and not the actual m/z values from the raw mass spectra.

(1035) Statistical significance of the results generated using the Q1-TOF and FBS-II MS was performed using the exact Cochran-Armitage test to find to compare the distributions of these specificity and sensitivity values between the two instrumental platforms evaluated since the models are constructed independently from each other.

Appendix A

Model 1	Test	Validity
Sensitivity	100%	100%
Specificity	100%	97%
Normal State	30/30 (100%)	38/37 (97%)
Ovarian Cancer State	57/57 (100%)	40/40 (100%)

Node	Count	State	StateSum	Entropy	m/z
0	7	1	7	0.000000	8688.074 8602.237 4644.793 7860.121 1464.593
1	3	0	0	0.000000	0 0.129292 1 0.404121 0.577349
2	6	1	6	0.000000	0 0.666673 1 0.230548 0.242727
3	16	1	16	0.000000	0 0.134574 1 0.891099 0.319833
4	3	0	0	0.000000	0 0.157213 1 0.091906 0.149974
5	13	1	11	0.000000	0 0.653332 0.714489 0.108038 1
6	4	0	0	0.000000	2 0.320183 1 0.123428 0.39002
7	2	1	2	0.000000	1 0.425972 1 0.178253 0.191267
8	2	1	2	0.000000	0 0.232833 1 0.146285 0.79188
9	2	0	0	0.000000	0 0.683164 0.613282 0.408828 1
10	2	1	2	0.000000	0 0.211845 0.666812 0.115333 1
11	5	0	0	0.000000	0 0.976017 0.954457 0.170028 0.628189
12	3	0	0	0.000000	1 0.341484 1 0.443244 0.367961
13	2	1	2	0.000000	0 0.14915 1 0.890447 0.340318
14	2	0	0	0.000000	0 0.682325 1 0.358043 0.558506
15	1	0	0	0.000000	0 0.859219 0.724638 0.26087 1
16	1	0	0	0.000000	0 0.645833 1 0.502063 0.835417
17	1	0	0	0.000000	0 0.794488 0.894737 0.894236 1
18	2	0	0	0.000000	0 0.97861 1 0.423498 0.63491
19	2	1	2	0.000000	0 0.448107 1 0.183052 0.753369

Model 2	Test	Validity
Sensitivity	100%	100%
Specificity	100%	85%
Normal State	30/30 (100%)	35/37 (95%)
Ovarian Cancer State	57/57 (100%)	40/40 (100%)

Node	Count	State	StateSum	Error	m/z
0	7	1	7	0	8605.678 5773.642 6258.91 7060.121 8708.065 748.048
1	3	0	0	0	0.836245 0.103495 0.112529 0.966826 0.445348 0
2	10	1	10	0	0.891918 0.304599 0.273147 0.468784 0.965088 0
3	3	0	0	0	0.069852 0.103221 0.545584 0.405998 0
4	13	1	8	0	0.668897 0.155636 0.241726 0.965208 0.984241 0
5	3	1	3	0	0.968501 0.107261 0.192038 0.625891 0.857142 0
6	2	0	0	0	0.585203 0.103657 0.125338 0 0.430678 0
7	3	1	3	0	0.810908 0.23603 0.555267 0.974007 0 0
8	8	1	8	0	0.294877 0.117567 0.231772 0 0.818855 0
9	7	0	0	0	0 0.112112 0.122808 0.745443 0.523186 0
10	10	1	10	0	0.59098 0.178288 0.258633 0.503651 0 0
11	1	0	0	0	0 0.047377 0.061628 0.284495 0.406995 0
12	4	0	0	0	0 0.133162 0.208333 0.305556 0.803241 0
13	1	1	1	0	0.59657 0.199346 0.30219 0.707978 0 0
14	1	0	0	0	0.0411785 0.12549 0.137255 0 0.266867 0
15	1	0	0	0	0.818951 0.311436 0.408759 0 0.981071 0
					0.885809 0.315909 0.404545 0.711364 0 0

Model 3	Test	Validity
Sensitivity	100%	100%
Specificity	100%	97%
Normal State	30/30 (100%)	36/37 (97%)
Ovarian Cancer State	57/57 (100%)	40/40 (100%)

Node	Count	State	StateSum	Error	m/z
0	30	1	29	1	8605.678 6668.674 6345.108 9843.343 3354.185 7374.687 5968.506
1	6	0	1	1	0.387228 0.015495 0.128241 0.008888 0.213168 0.032584
2	5	0	0	1	0.813858 0.049105 0.268494 0.031082 0.354791 0.060409
3	18	1	14	5	0.943076 0.9957 0.023126 0.32079 0.05742 0.600263 0.033526
4	1	0	0	0	0.582078 0.049422 0.20029 0.026914 0.389413 0.028103
5	1	0	0	0	0.918869 0.042514 0.260628 0.170055 0.914972 0
6	3	1	3	0	0.620513 0.125356 0 0.333333 0.848718 0.321937
7	1	1	1	0	0.715204 0.006153 0.19088 0.080695 0.722323 0.025888
8	3	0	0	0	0.573182 0.0203 0.151675 0.130511 0.082363 0.044092
9	3	0	0	0	0.937262 0.9936 0.115137 0.159158 0 0.830834 0.113328
10	11	0	0	0	0.0722109 0.017883 0.045724 0.057432 0.817682 0.059098
11	2	1	2	0	0.850943 0.320755 0.230189 0 0.684151 0.301687
12	1	0	0	0	0.041404 0.078637 0.148901 0.038538 0.845357 0
13	1	0	0	0	0.980798 0.075332 0.51551 0.0401773 0.025111 0
					0.908907 0.081081 0.012012 0.189189 0.429429 0

Model 4	Test	Validity
Sensitivity	100%	100%
Specificity	100%	100%
Normal State	30/30 (100%)	37/37 (100%)
Ovarian Cancer State	57/57 (100%)	40/40 (100%)

Node	Count	State	StateSum	Error
0	8	1	8	1
1	3	0	0	0
2	10	1	10	1
3	3	0	0	0
4	8	1	8	1
5	10	1	10	1
6	3	0	0	0
7	4	1	4	1
8	2	0	0	0
9	1	1	1	1
10	8	0	0	0
11	10	1	10	1
12	3	0	0	0
13	3	0	0	0
14	2	0	0	0
15	1	1	1	1

m/z	7060.121	7096.622	8605.678	8548.771	8708.085	818.4801	8540.636	8352.723
0	0.917113	0.21551	0.861398	0.121208	0.444445	0	0.518113	0.110812
0	0.492091	0.305348	0.866398	0.205158	0.894171	0	0.951383	0.236869
0	0.547569	0.173663	1	0.104231	0.409816	0	0.51695	0.092858
0	0.929844	0.33378	0.874226	0.166695	0.863815	0	0.90104	0.157423
0	0.732832	0.276296	1	0.135825	0.570398	0	0.883495	0.107333
3	0.846923	0.304081	0.883209	0.148318	0.82462	0	0.918506	0.12435
0	0.346591	0.221128	1	0.173951	0.806024	0	0.827509	0.179187
0	1	0.262028	0.56594	0.124258	0.40729	0	0.422331	0.10647
0	0.794377	0.531631	0.515963	0.290957	0.814304	0	1	0.29789
0	1	0.270158	0.932108	0.145688	0.831683	0	0.946252	0.132958
0	0.437313	0.281307	0.815518	0.170126	0.890082	0	0.886262	0.143119
0	0.282368	0.113517	1	0.06052	0.405555	0	0.507878	0.047164
0	0.652298	0.545467	0.758154	0.391447	0.893289	0	0.878634	0.381204
0	0.663094	0.35973	0.601834	0.214181	0.872875	0	1	0.191813
1	1	0.638476	0.845795	0.372277	0.857743	0	0.965217	0.311208
0	1	0.237164	0.735178	0.105402	0.753623	0	0.756258	0.102767

W01005011474

PCTUS00402413

Model 5	Test	Validity
Sensitivity	100%	100%
Specificity	100%	97%
Normal State	30/30 (100%)	36/37 (97%)
Ovarian Cancer State	57/57 (100%)	40/40 (100%)

Node	Count	State	StateSum	Error
0	30	1	30	1
1	2	0	0	0
2	2	0	0	0
3	17	1	17	1
4	2	0	0	0
5	5	1	5	1
6	1	0	0	0
7	2	1	2	1
8	2	0	0	0
9	1	0	0	0
10	2	0	0	0
11	3	0	0	0
12	2	0	0	0
13	1	0	0	0
14	1	0	0	0
15	1	0	0	0
16	1	0	0	0
17	2	1	2	1

m/z	11601.83	8718.517	3449.205	4260.403	1229.752	2007.145	8602.237	7060.121	848.10
0	0.045973	0.188825	0.031338	0.084657	0.008804	0.010191	1	0.232181	0.0142
0	0.190458	0.752349	0.206444	0.438551	0	0.0639	1	0.321633	0.3785
0	0.185637	0.728544	0.16697	0.355362	0	0.029894	0.730038	1	0.052
6	0.078896	0.33797	0.088988	0.20709	0.028185	0.022459	1	0.437262	0.0432
0	0.115091	0.512947	0.110247	0.353818	0.002048	0.043823	1	0.230498	0.2096
0	0.090591	0.287811	0.087215	0.154745	0.018448	0.049325	1	0.740332	0.0142
0	0.202229	0.542884	0.402866	0.52707	0.197452	0	0.821019	1	0.258
0	0.108417	0.226812	0.165819	0.205581	0.014039	0.018811	0.69384	1	0.035
0	0.143113	1	0.214746	0.826275	0.086988	0	0.82163	0.582268	0.483
0	0.178571	0.921053	0.274438	0.744361	0	0.087669	1	0.772556	0.24
0	0.127322	0.855385	0.288389	0.341074	0.000843	0.066154	0.973585	0.601901	0.558
0	0.230129	0.728008	0.280667	0.633693	0.045805	0.024148	0.754434	1	0.104
0	0.18007	0.762553	0.208338	0.57439	0	0.086841	1	0.875463	0.400
0	0.127701	0.665815	0.125737	0.875835	0.037328	0	1	0.844794	0.149
0	0.138095	0.784127	0.163492	0.477778	0	0.014288	1	0.760317	0.083
0	0.291045	0.808458	0.271144	0.41791	0	0.014825	0.895522	1	0.383
0	0.158163	0.785714	0.318878	0.558673	0	0.035714	1	0.812245	0.277
0	0.154471	0.472129	0.131158	0.216488	0.027697	0	1	0.784209	0.187

W01005011474

PCTUS00402413

Model 8	Test	Validity
Sensitivity	100%	100%
Specificity	100%	97%
Normal State	30/30 (100%)	36/37 (97%)
Ovarian Cancer State	57/57 (100%)	40/40 (100%)

Node	Count	State	StateSum	Error	m/z
0	12	1	12	0	8688.674 8602.237 7060.121 4920.131 10431.02 2817.487
1	2	0	0	0	0 0.212058 1 0.44328 0.05893 0.243358 0
2	19	1	19	0	0 0.71195 1 0.320383 0.194065 0.325602 0
3	6	0	0	0	0 0.181351 1 0.188047 0.02488 0.074401 0
4	7	1	5	2	0 0.721687 0.728508 1 0.146458 0.244383 0
5	9	1	6	3	0 0.328961 1 0.382833 0.054395 0.118492 0
6	4	0	0	0	0 0.430797 1 0.448852 0.061423 0.253657 0
7	3	1	3	0	0 0.479383 1 0.241389 0.13775 0.184372 0
8	1	1	1	0	0 0.265618 1 0.781812 0.070789 0.199972 0
9	1	1	1	0	0 0.264708 0.703013 1 0.066715 0.351509 0
10	6	0	0	0	0 0.218579 1 0.872131 0.213115 0.464481 0
11	2	0	0	0	0 0.979239 0.960158 0.898689 0.134247 0.169243 0
12	1	0	0	0	0 0.687882 1 0.587495 0.248281 0.240037 0
13	1	0	0	0	0 0.195426 0.60499 1 0.04262 0.096074 0
14	1	0	0	0	0 0.686347 1 0.254244 0.158827 0.560888 0
15	1	0	0	0	0 0.788458 0.890625 1 0.330729 0.58625 0
16	1	0	0	0	0 0.987805 1 0.536585 0.140244 0 0
17	1	1	1	0	0 0.486765 1 0.741176 0.068177 0.448529 0
					0 0.478368 1 0.886279 0.089999 0.25859 0

W0.J05011174

PCTUS000102413

Model 8	Test	Validity
Sensitivity	100%	100%
Specificity	100%	95%
Normal State	30/30 (100%)	35/37 (95%)
Ovarian Cancer State	57/57 (100%)	40/40 (100%)

Node	Count	State	StateSum	Error	m/z
0	9	1	9	0	708.1657 8605.678 6806.648 7060.121 6761.677 2472.108 8708.095 5511.917 1195.325 500
1	15	0	0	0	0 0.978759 0.129335 0.890028 0.141874 0.08436 0.465115 0.117084 0.112831 0.01
2	15	1	15	0	0 0.994064 0.168514 0.384269 0.247893 0.076075 0.898872 0.147354 0.126049 0.11
3	4	0	0	0	0 0.092694 0.597216 0.154853 0.061148 0.463791 0.081717 0.104318 0.03
4	12	1	8	4	0 0.660345 0.19312 0.987633 0.301109 0.102143 0.97033 0.184698 0.154734 0.11
5	4	1	4	0	0.008896 0.966228 0.160728 0.835588 0.230458 0.048255 0.860368 0.09372 0.147285 0.0
6	1	0	0	0	0 0.548765 0.094072 1 0.130738 0.048314 0.384022 0.087314 0.084237 0.0
7	1	1	1	0	0 0.589539 0.283537 0.972581 0.705793 0.10061 1 0.181402 0.385871 0.2
8	3	1	3	0	0 0.807692 0.048154 1 0.084815 0.181538 0.423077 0.038462 0.315385 0.2
9	5	0	0	0	0 0.882666 0.160095 1 0.274763 0.083765 0.814652 0.091038 0.151456 0.1
10	10	1	10	0	0 0.67702 0.16947 0.448973 0.283484 0.093472 1 0.116756 0.184678 0.1
11	2	0	0	0	0.001145 1 0.082802 0.272652 0.076581 0.027031 0.397883 0.035259 0.049178 0.0
12	4	0	0	0	0 0.701871 0.326662 0.583559 0.401201 0.083416 1 0.270312 0.134082 0.2
13	1	0	0	0	0 0.585976 0.201684 0.698887 0.327029 0.059885 1 0.153016 0.12643 0.1
14	1	0	0	0	0 0.810258 0.305128 1 0.412821 0.002584 0.958974 0.269231 0.010256 0.2
					0 0.8742 0.347548 0.729211 0.663113 0.132186 1 0.289979 0.249487 0.2

W0.J05011173

PCTUS000102413

Model 9	Test	Validity
Sensitivity	100%	100%
Specificity	100%	100%
Normal State	30/30 (100%)	37/37 (100%)
Ovarian Cancer State	57/57 (100%)	40/40 (100%)

Node	Count	State	StateSum	Error	m/z
0	28	1	28	0	7048.018 8602.237 8684.385 1144.796 4260.403
1	4	0	0	0	0 0.117795 0.189138 0.00018 0.098848
2	3	0	0	0	0 0.44898 0.724911 0 0.518048
3	12	1	9	0	0 0.618286 0.893434 0.814925 0 0.472577
4	7	0	1	0	0 0.191145 0.325061 0 0.163693
5	9	1	9	0	0 0.214739 0.50704 0 0.340581
6	4	0	0	0	0 0.3496 0.388951 0 0.221401
7	1	0	0	0	0 0.745345 0.698562 0 0.634987
8	1	1	1	0	0 0.740741 0.618519 0 0.522222
9	1	0	0	0	0 0.846484 0.373047 0 0.303711
10	2	0	0	0	0 0.48337 0.846888 0 0.897438
11	1	0	0	0	0 0.515608 0.903218 0 0.728896
12	1	0	0	0	0 0.739758 0.882573 0 0.944444
13	1	1	1	0	0 0.512546 0.25980 0 0.108527
14	1	0	0	0	0 0.348457 0.802382 0 0.675197
15	1	0	0	0	0 0.833148 0.793872 0 0.485181

Model 10	Test	Validity
Sensitivity	100%	100%
Specificity	100%	100%
Normal State	30/30 (100%)	37/37 (100%)
Ovarian Cancer State	57/57 (100%)	40/40 (100%)

Node	Count	State	StateSum	Error	m/z
0	4	1	4	0	7202.718 6004.417 7080.121 1001.654 1255.593 8587.113 4377.854 8605.678 8709.5
1	2	0	0	0	0 0.173188 0.074963 0.970482 0.003208 0.042688 0.283611 0.14722 0.858894 0.371
2	10	1	10	0	0 0.319725 0.176894 0.393018 0 0.164671 0.825969 0.378272 0.917131 0.8672
3	2	0	1	0	0 0.189442 0.082052 0.680658 0 0.055131 0.403148 0.151314 0 0.4594
4	2	1	2	0	0 0.361857 0.113665 0 0 0.121266 0.562191 0.202878 0.70216 0.8284
5	1	1	1	0	0 0.213108 0.072628 0.578897 0 0.050348 0.662743 0.155164 0 0.5021
6	3	1	3	0	0 0.284091 0.113638 0.940341 0 0.150588 0.805114 0.207386 0 0.4711
7	7	1	5	0	0 0.263962 0.121837 0.831316 0 0.080509 0.411379 0.183044 0 0.6011
8	2	1	2	0	0 0.235242 0.08713 0.676821 0 0.082517 0.506815 0.140705 0 0.8661
9	2	1	2	0	0 0.081188 0.421919 0.159605 0.819174 0.3851
10	1	0	0	0	0 0.227143 0.128887 0 0 0.068565 0.418376 0.128141 0.52401
11	1	0	0	0	0 0.280298 0.087375 0.746658 0 0.15758 0.302414 0.123253 0.472681
12	1	0	0	0	0 0.584188 0.180432 0.791814 0 0.174551 0.597064 0.17292 0.882055
13	1	1	1	0	0 0.383381 0.168028 0.71915 0 0.04466 0.188106 0.105068 0.483184 0.4301
14	2	1	2	0	0 0.254143 0.094835 0 0 0.088878 0.388488 0.190463 0 0.822
15	2	1	2	0	0 0.484788 0.101004 0.847496 0 0.083505 0.313402 0.130928 0 0.904
16	1	1	1	0	0 0.303093 0.053808 0.485979 0 0 0.125874 0.454545 0.202787 0.825175 0.573
17	1	1	1	0	0 0.237782 0.167832 0 0 0 0.070398 0.522135 0.262555 0.833444 0.971
18	2	0	0	0	0 0.335049 0.16409 0.489544 0 0 0.105538 0.508054 0.173701 0.830654 0.874
19	2	0	0	0	0 0.369559 0.068285 0 0 0 0.106513 0.341438 0.109485 0.518447
20	2	0	0	0	0 0.243242 0.067837 0.335432 0 0 0.045892 0.286063 0.113572 0 0.382
21	8	1	8	0	0 0.123575 0.048128 0.211115 0 0 0.113593 0.450127 0.132826 0.790771
22	2	0	0	0	0 0.211598 0.059312 0.548008 0 0 0 0 0 0

21	4	0	0	0	0.329776	0.110944	0.506651	0	0.132027	0.484959	0.18387	0.587533
22	1	0	0	0	0.253837	0.126328	0.291617	0	0.11098	0.5183	0.26307	0.9185
23	1	0	0	0	0.601351	0.344595	0.783514	0	0.098847	0.88036	0.481882	0.878378
24	1	0	0	0	0.329101	0.118402	0.569312	0	0.078181	0.274074	0.111111	0.394709
25	2	0	0	0	0.453461	0.170665	0.800839	0	0.118823	0.618038	0.254696	0.552077
26	3	1	3	0	0.119065	0.10091	0.491402	0	0.082836	0.204372	0.145723	0.2959
27	1	0	0	0	0.178475	0.119283	0.300448	0	0.101345	0.817489	0.220828	0.873543
28	1	0	0	0	0.554658	0.297571	0.870445	0	0.109312	0.534413	0.317814	0.720848
29	1	1	1	0	0.083584	0.030732	0.097721	0	0.02787	0.11882	0.058356	0.3080
30	1	0	0	0	0.457023	0.180284	0.57652	0	0.125788	0.574423	0.400419	0.698113
31	1	0	0	0	0.678325	0.278371	0.796287	0	0.187784	0.601256	0.398734	0.879747
32	1	1	1	0	0.189982	0.060578	0.289331	0	0.063291	0.352622	0.136528	0.6084

WO 2005011473

PGT0500402410

Model 11	Test	Validity
Sensitivity	100%	100%
Specificity	100%	97%
Normal State	30/30 (100%)	95/97 (97%)
Ovarian Cancer State	57/57 (100%)	40/40 (100%)

Node	Count	State	StateSum	Error	m/c
0	5	1	5	0	0.443439
1	1	0	0	0	0.005091
2	2	1	2	0	0.129668
3	1	1	1	0	0.02943
4	3	0	0	0	0.323752
5	6	1	6	0	0.192401
6	1	1	1	0	0.184719
7	2	1	2	0	0.12639
8	4	1	4	0	0.12784
9	3	1	3	0	0.11335
10	1	0	0	0	0.120282
11	1	1	1	0	0.124895
12	2	0	0	0	0.10223
13	1	1	1	0	0.106281
14	4	1	4	0	0.14548
15	1	1	1	0	0.153571
16	2	1	2	0	0.139768
17	2	1	2	0	0.157544
18	1	1	1	0	0.184549
19	1	0	0	0	0.128781
20	1	1	1	0	0.150685

WO 2005011473

PGT0500402410

21	2	0	0	0	0.0325044	0.8441	0.868557	0	1	0.337625	0.940181	0.381157	0.190878	0.33778
22	7	1	6	1	0.182207	1	0.306461	0	0.651426	0.178093	0.223762	0.212152	0.060816	0.19411
23	3	0	0	0	0.185007	1	0.573712	0	0.801758	0.254952	0.718139	0.318285	0.190385	0.33331
24	1	0	0	0	0.0512298	1	0.663598	0	0.855263	0.532895	0.434211	0.475877	0.236842	0.34421
25	2	1	2	0	0.154238	1	0.55516	0.033829	0.433112	0.279879	0.2783	0.612885	0.009195	0.2719
26	2	0	0	0	0.023325	0.812591	0.702928	0	1	0.325657	0.802751	0.403117	0.171496	0.32541
27	1	0	0	0	0.2737	1	0.227829	0	0.840673	0.30581	0.440367	0.366972	0.114579	0.24811
28	1	1	1	0	0.244582	1	0.547988	0	0.668731	0.390093	0.340557	0.804954	0.065018	0.31261
29	1	0	0	0	0.357798	1	0.443425	0	1	0.357798	0.82263	0.477084	0.140673	0.339
30	1	0	0	0	0.00255	0.897996	0.484481	0	1	0.48818	0.741348	0.890893	0.127505	0.33871
31	2	0	0	0	0.285715	0.975111	0.493958	0	0.856469	0.375087	0.519177	0.279772	0.161407	0.39231
32	1	0	0	0	0.225316	0.920886	0.300833	0	0.893671	0.389241	0.512858	1	0.018987	0.51891
33	3	1	3	0	0.087097	1	0.183182	0	0.445501	0.125678	0.182821	0.453178	0.055861	0.1002
34	1	1	1	0	0.078783	1	0.144424	0	0.26891	0.086892	0.122486	0.025584	0.074954	0.0678
35	1	0	0	0	0.02891	0.884928	0.846919	0	1	0.400474	0.599526	0.575829	0.151659	0.2725
36	1	1	1	0	0.178923	1	0.517028	0	0.377709	0.227554	0.308598	0.28483	0.085139	0.3287
37	1	1	1	0	0.030443	1	0.612768	0	0.859574	0.282979	0.353191	0.653182	0.131915	0.3914
38	1	0	0	0	0.030571	0.748018	0.565478	0	1	0.328365	0.414683	0.337302	0.188492	0.1984

W/O 2005011474

PCT/US2004/024413

Model 12	Test	Validity
Sensitivity	100%	100%
Specificity	100%	95%
Normal State	30/30 (100%)	35/37 (95%)
Ovarian Cancer State	67/67 (100%)	40/40 (100%)

Node	Count	State	StateSum	End	m/z
0	6	1	6	1	8885.2 : 8709.548 : 7065.771 : 1132.049 : 8603.678
1	2	0	1	1	0.227356 : 0.285099 : 0.294878
2	5	1	5	1	1.0579418 : 0.996678 : 0.249831
3	2	0	0	0	0.288212 : 0.48104 : 0.337354
4	2	1	2	1	0.839955 : 1 : 0.945907
5	7	1	7	1	0.444594 : 0.494724 : 0.255931
6	8	1	8	1	0.328118 : 0.404857 : 0.471829
7	6	1	6	1	0.420976 : 0.599319 : 0.470768
8	3	0	0	0	0.81664 : 0.802203 : 0.35835
9	1	1	1	1	0.653035 : 0.84379 : 0.223522
10	4	0	0	0	0.645 : 0.645 : 0.9678
11	1	0	0	0	0.430854 : 1 : 0.405585
12	11	1	11	1	0.155009 : 1 : 0.448905
13	1	1	1	1	0.281647 : 0.357538 : 0.14863
14	1	1	1	1	0.650505 : 1 : 0.39596
15	1	1	1	1	0.313343 : 0.812594 : 1
16	2	1	2	1	0.640593 : 0.804083 : 0.442778
17	1	0	0	0	0.771379 : 1 : 0.319372
18	2	1	2	1	0.395313 : 0.746381 : 0.548285
19	2	0	0	0	0.358261 : 1 : 0.141059
20	2	0	0	0	0.357038 : 1 : 0.251898
	1	0	0	0	0.986008 : 1 : 0.68272
					0 : 0.847026

W/O 2005011474

PCT/US2004/024413

21	1	0	0	0.0334625	1	0.31137	0	0.260982	
22	1	1	1	0.0376208	0.533762	1	0	0.951768	
23	2	0	0	0.0358085	1	0.272823	0	0.537859	
24	2	0	0	0.0579131	1	0.240333	0	0.640437	
25	1	0	0	0.0471058	1	0.660679	0	0.51497	
26	1	0	0	0.066581	1	0.398458	0	0.62982	
27	1	1	1	0.019256	0.833698	0.669584	0	1	
28	1	0	0	0.0782258	1	0.628032	0	0.846774	
29	1	1	1	0	0.516	1	0.518	0	0.896
30	1	1	1	0.0403558	0.594569	0.152622	0	1	

WO 2005011474

PC/US2004/074413

Model:13	Test:	Validity:
Sensitivity	100%	100%
Specificity	100%	95%
Normal State	30/30 (100%)	35/37 (95%)
Ovarian Cancer State	57/57 (100%)	40/40 (100%)

WO 2005011474

Node	Count	State	StateSum	Error	m/z
0	0	1	0	0	11098.07 6501.789 2087.371 8805.678 8688.674 7048.838 4262.107
1	1	0	0	0	0.053842 0.050306 0 1 0.277113 0.258017 0.126978
2	1	0	0	0	0.194368 0.016801 0 1 0.780282 0.24507 0.416901
3	1	0	0	0	0.230024 0.178177 0 1 0.890315 0.738077 0.483947
4	1	0	0	0	2.047783 -0.03069 0.000757 1 0.473931 0.24506 0.11983
5	1	0	0	0	0.074636 0.064462 0 1 0.43221 0.343755 0.20137
6	1	0	0	0	0.094925 0.130789 0 1 0.671894 0.378017 0.273367
7	1	0	0	0	0.059567 0.032491 0 1 0.644404 0.355596 0.034296
8	1	0	0	0	0.230797 0.139693 0 1 0.830324 0.189319 0.459966
9	1	0	0	0	0.205333 0.056 0 1 0.514687 0.794687 0.122867
10	1	0	0	0	0.108929 0.123214 0 0 0.0921429 1 0.883829 0.457143
11	1	0	0	0	0.068063 0.108377 0 0.832461 0.997382 1 0.505238
12	1	0	0	0	0.0378 0.018129 0.005735 1 0.292722 0.108974 0.075537
13	1	0	0	0	0.068486 0.115332 0 0.82768 0.49322 1 0.238806
14	1	0	0	0	0.082474 0.195878 1 0.402082 0.237113 0.154839
15	1	0	0	0	0.12328 0.280316 0 1 0.852883 0.274354 0.310139
16	1	0	0	0	0.043462 0.066573 0 1 0.835869 0.380821 0.514702
17	1	0	0	0	0.124457 0.059334 0 1 0.609262 0.357453 0.444284
18	1	0	0	0	0.192394 0.127517 0 0.876957 1 0.438479 0.628635
19	1	0	0	0	0.091245 0.165228 0 1 0.641184 0.181258 0.282387
20	1	1	1	0	0.0313726 0.124183 0.95098 1 0.650327 0.441178
					0.153302 0.178245 0 1 0.415094 0.566038 0.235849

PC/US2004/074413

21	1	0	0	0	0.128719	0.185842	0	0.759901	1	0.675743	0.537129
22	2	0	0	0	0.194312	0.20655	0	0.84264	1	0.528225	0.430212
23	1	1	1	0	0.2125	0.2	0	0.905	0.47	1	0.19
24	1	0	0	0	0.270089	0.084821	0	0.841518	1	0.870536	0.546875
25	1	0	0	0	0.134441	0.128399	0	0.880383	1	0.311178	0.303825
26	1	0	0	0	0.387436	0.339744	0	0.858974	1	0.903846	0.490385
27	1	0	0	0	0	0.257908	0	0.824574	1	0.491484	0.593874
28	1	0	0	0	0.29085	0.382745	0	1	0.973856	0.990195	0.470588
29	1	0	0	0	0	0.147287	0.038176	0.976744	1	0.50846	0.423773
30	1	0	0	0	0.047222	0.175	0	0.76	1	0.497222	0.480556
31	1	1	1	1	0.18896	0.278658	0	1	0.733202	0.743063	0.320158
32	1	0	0	0	0.061404	0.285088	0	0.813596	1	0.598684	0.33114
33	1	1	1	1	0.090909	0.130165	0	1	0.733471	0.607438	0.208678

WO 2005011374

PCT/US2004/024413

Model 14	Test	Validity
Sensitivity	100%	100%
Specificity	100%	97%
Normal State	30/30 (100%)	36/37 (97%)
Ovarian Cancer State	57/57 (100%)	40/40 (100%)

Node	Count	State	StateSum	Entropy	m/z
0	3	1	3	1	4162.719 8588.487 8708.548 8664.385 1318.958 8805.678 2280.258 7060.121
1	1	0	0	0	0.095692 0.344856 0.318228 0.242556 0.007524 0.969059 0.000948 0.959932
2	5	1	5	0	0.488175 0.688894 1 0.626728 0 0.880184 0.004608 0.31106
3	6	1	6	0	0.117272 0.439504 0.401233 0.30528 0 1 0.039692 0.653983
4	1	0	0	0	0.083016 0.498557 0.325581 0.28407 0.00115 1 0.014817 0.410254
5	1	1	1	0	0.153971 0.58671 0.85624 0.684506 0 0.662885 0.008483 1
6	3	1	3	0	0.108524 0.591887 0.504762 0.657143 0 1 0.105952 0.55
7	2	1	2	0	0.127988 0.493341 0.417544 0.3649 0.002772 0.984158 0.050381 0.925263
8	7	1	7	0	0.207404 0.724887 0.602076 0.532475 0 1 0.037808 0.814917
9	1	0	0	0	2 0.178899 0.716138 0.812647 0.551972 0.005477 0.998362 0.018468 0.650556
10	1	1	1	0	0.687282 0.824477 0.827697 0.65599 0 1 0.119163 0.310789
11	2	1	2	0	0.108787 0.426778 0.381227 0.403068 0 0.559275 0.026499 1
12	3	0	0	0	0.108972 0.628005 0.453237 0.653568 0.005034 1 0.030471 0.406813
13	1	1	1	0	0.152024 0.439361 1 0.428457 0.005728 0.478396 0.0065 0.730046
14	2	1	2	0	0.109208 0.304069 0.432548 0.248253 0 0.441114 0.068523 1
15	1	1	1	0	0.253559 0.857705 0.691482 0.692784 0.013308 1 0.006591 0.448839
16	1	1	1	0	0.242188 0.335938 0.623438 0.826125 0 0.804688 0.226562 1
17	1	0	0	0	0.225275 0.807892 1 0.723443 0.021978 0.908425 0 0.448718
18	1	1	1	0	0.162909 0.687706 0.880555 0.605697 0.043478 0.828036 0 1
19	2	0	0	0	0.142689 0.621053 0.768421 0.492398 0.014035 1 0 0.817544
20	5	1	5	0	0.172991 0.469998 1 0.484749 0.004408 0.484017 0 0.287822
					0.082151 0.474033 0.407828 0.324867 0 1 0.013184 0.257672

WO 2005011374

PCT/US2004/024413

21	2	0	0	0	0.16018	0.506442	1	0.439991	0.008219	0.7738	0	0.511529
22	3	1	3	0	0.153858	0.858383	0.460659	0.432758	0.004074	1	0.033124	0.717648
23	1	1	1	0	0.2021	0.645689	0.703412	0.671916	0.026247	1	0	0.55843
24	4	0	0	0	0.2007	0.575851	1	0.530549	0	0.522831	0.024103	0.458878
25	1	0	0	0	0.209799	0.757538	0.913317	0.504271	0	1	0.035178	0.246231
26	2	0	0	0	0.387106	0.8472	1	0.835188	0	0.850562	0.070583	0.702618
27	1	1	1	0	0.164818	0.438886	0.28794	0.282092	0	0.728002	0.041204	1
28	1	0	0	0	0.132353	0.438914	1	0.535973	0	0.352941	0.001131	0.539893
29	1	1	1	0	0.123829	0.300728	0.240375	0.207078	0.009385	0.37593	0	1
30	2	0	0	0	0.222129	0.625428	1	0.575785	0	0.504059	0.049331	0.778349
31	1	0	0	0	0.101895	0.62343	1	0.57328	0	0.637089	0.041874	0.222333
32	1	0	0	0	0.232258	0.673118	1	0.812903	0	0.703226	0.124731	0.862366
33	2	1	2	0	0.132722	0.535895	0.83435	0.388105	0.008025	1	0.01613	0.559726
34	1	1	1	0	0.035639	0.539873	0.292872	0.246295	0	1	0.000708	0.077982
35	1	0	0	0	0.306122	0.716837	1	0.665816	0	0.632853	0.030612	0.484894
36	1	1	1	0	0.210428	0.724395	0.787709	0.581006	0	0.929236	0.130354	1
37	1	1	1	0	0.154391	0.627479	0.787835	0.405069	0.031162	0.715297	0	1
38	1	1	1	0	0.070746	0.626195	0.586042	0.378585	0	1	0.016296	0.248566

W/O 2005011674

PCVUS004/024113

Model 15	Test	Validity
Sensitivity	100%	100%
Specificity	100%	100%
Normal State	30/30 (100%)	37/37 (100%)
Ovarian Cancer State	57/57 (100%)	40/40 (100%)

Node	Count	State	StateSum	Error	m/z
0	33	1	33	1	8870.938 2374.244 1278.881 7060.121 4292.9 8706.065 8605.678
1	23	1	16	7	0.120039 0.024623 0.011125 0.949945 0.171834 0.527519 0.872824
2	7	0	2	5	0.141853 0.02381 0.020886 0.628684 0.162889 0.626018 0.899723
3	16	0	1	15	0.186489 0 0.153321 0.882675 0.152271 0.953348 0.714632
4	3	1	3	0	0.144659 0 0.181107 0.595845 0.178005 1 0.741838
5	1	1	1	0	0.056997 0 0.043224 1 0.088753 0.359943 0.468551
6	1	0	0	1	0.04065 0 0.000353 0.076362 0.138211 0.276423 1
					0.358639 0.146597 0 0.337696 0.397806 1 0.984293

W/O 2005011674

PCVUS004/024113

What is claimed is:

1. A model usable in determining whether a biological sample taken from a subject indicates that the subject has ovarian cancer, comprising:
a vector space having at least three dimensions; and
at least one diagnostic cluster defined in said vector space; said diagnostic cluster corresponding to one of a diseased cluster and a healthy cluster;
said vector space having a first dimension that corresponds to a first mass to charge ratio value from a mass spectrum; said first mass to charge ratio being about 7060;
said vector space having a second dimension that corresponds to a second mass to charge ratio value from a mass spectrum; said second mass to charge ratio being about 8605; and
said vector space having a third dimension that corresponds to a third mass to charge ratio value from a mass spectrum; said third mass to charge ratio being about 8706.
2. The model of claim 1, wherein the vector space has at least four dimensions; said vector space having a fourth dimension that corresponds to a fourth mass to charge ratio value from a mass spectrum; said fourth mass to charge ratio being about 5548.
3. A model usable in determining whether a biological sample taken from a subject indicates that the subject has ovarian cancer, comprising:
~~a vector space having at least three dimensions; and~~
at least one diagnostic cluster defined in said vector space; said diagnostic cluster corresponding to one of a diseased cluster and a healthy cluster;
said vector space having a first dimension that corresponds to a first mass to charge ratio value from a mass spectrum; said first mass to charge ratio being about 9807;
said vector space having a second dimension that corresponds to a second mass to charge ratio value from a mass spectrum; said second mass to charge ratio being about 2374; and
said vector space having a third dimension that corresponds to a third mass to charge ratio value from a mass spectrum; said third mass to charge ratio being about 1276.

4. The model of claim 3, wherein the vector space has at least four dimensions; said vector space having a fourth dimension that corresponds to a fourth mass to charge ratio value from a mass spectrum; said fourth mass to charge ratio being about 4292.
5. A method of determining whether a biological sample taken from a subject indicates that the subject has ovarian cancer by analyzing the biological sample to obtain a data stream that describes the biological sample, comprising:
a. abstracting the data stream to produce a sample vector that characterizes the data stream in a predetermined vector space containing a diagnostic cluster; the diagnostic cluster being an ovarian cancer cluster; the ovarian cancer cluster corresponding to the presence of ovarian cancer;
b. determining whether the sample vector rests within the ovarian cancer cluster; and
c. if the sample vector rests within the ovarian cancer cluster, identifying the biological sample as being taken from a subject that has ovarian cancer.

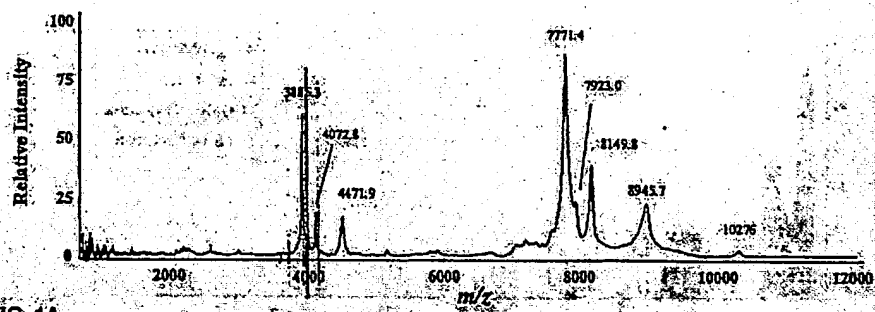


FIG. 1A

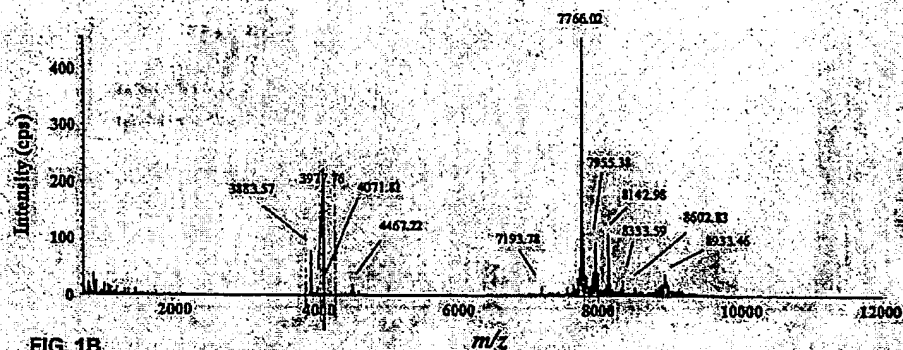


FIG. 1B

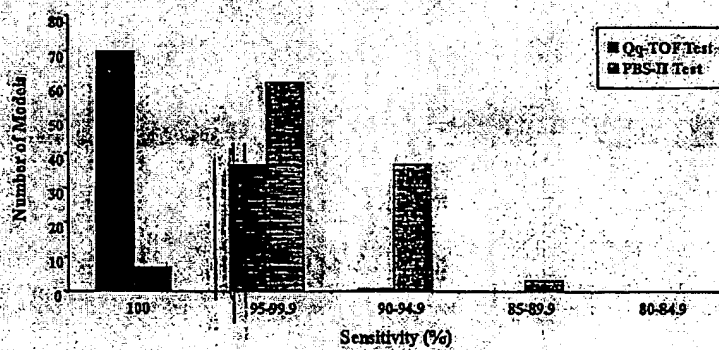


FIG. 2A

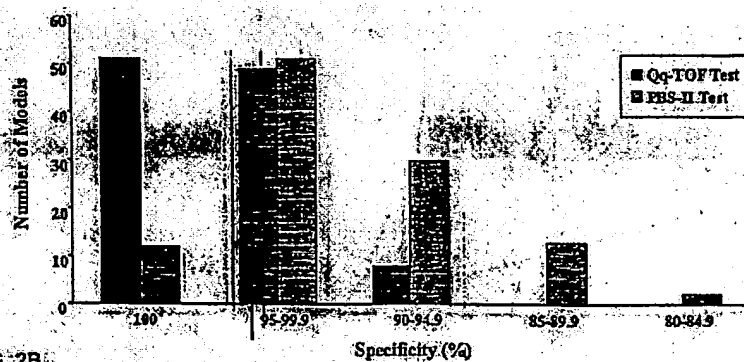


FIG. 2B

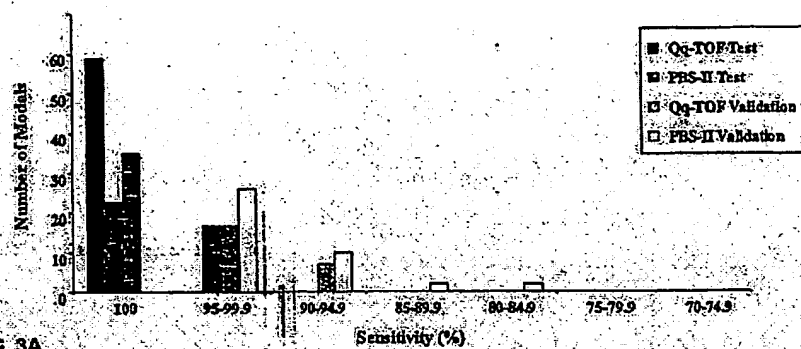


FIG. 3A

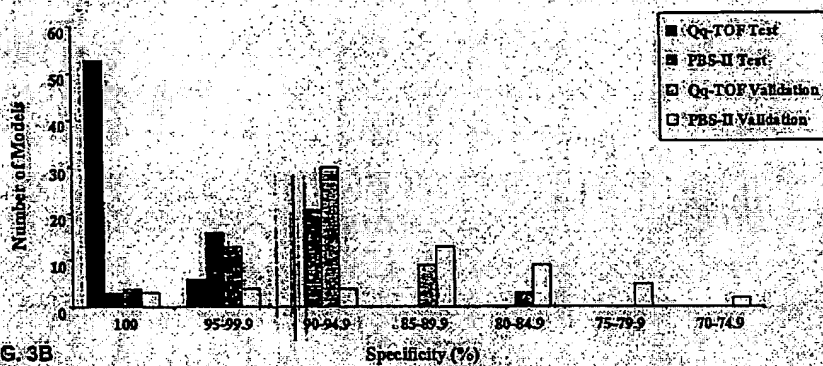


FIG. 3B

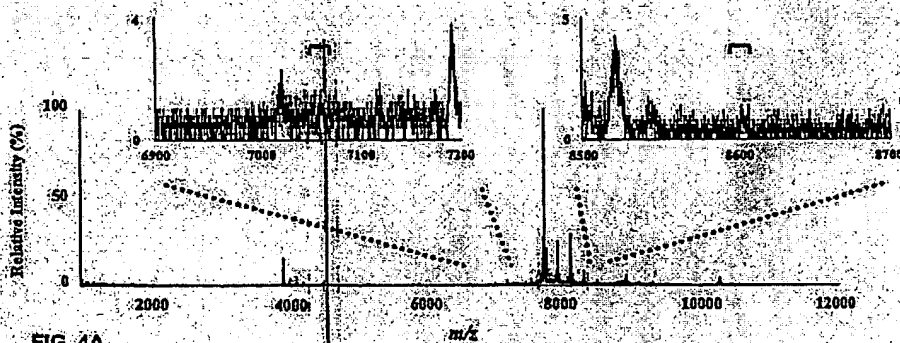


FIG. 4A

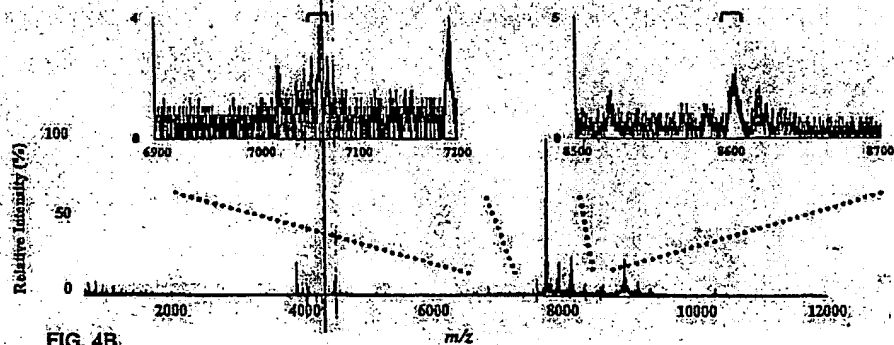


FIG. 4B

WO 2006/011413

PCT/US2006/011413

WO 2006/011413

PCT/US2006/011413

**This Page is Inserted by IFW Indexing and Scanning
Operations and is not part of the Official Record**

BEST AVAILABLE IMAGES

Defective images within this document are accurate representations of the original documents submitted by the applicant.

Defects in the images include but are not limited to the items checked:

- ☐ **BLACK BORDERS**
- ☐ **IMAGE CUT OFF AT TOP, BOTTOM OR SIDES**
- ☐ **FADED TEXT OR DRAWING**
- ☐ **BLURRED OR ILLEGIBLE TEXT OR DRAWING**
- ☐ **SKEWED/SLANTED IMAGES**
- ☐ **COLOR OR BLACK AND WHITE PHOTOGRAPHS**
- ☐ **GRAY SCALE DOCUMENTS**
- ☐ **LINES OR MARKS ON ORIGINAL DOCUMENT**
- ☐ **REFERENCE(S) OR EXHIBIT(S) SUBMITTED ARE POOR QUALITY**
- ☐ **OTHER:** _____

IMAGES ARE BEST AVAILABLE COPY.

As rescanning these documents will not correct the image problems checked, please do not report these problems to the IFW Image Problem Mailbox.